

24-Hour Pediatric Elbow

Ergonomic Control Button
Takes only single move to manipulate the suction source.

Soft Sleeve
Fully covers the suction catheter for better protection on both healthcare practitioners and the patients. Meanwhile, the transparency and the softness of the material offer easy operation for the procedure.

One-way Irrigation Port
Helps not only maintain the pressure of the ventilator circuit but also takes the injection of normal saline to clean up the catheter after each suction procedure.


BIOTEQ® 24-Hour Pediatric Elbow Catheter set is designed for safely suction neonatal and pediatric patients on mechanical ventilation, elbow adapter for standard and high-frequency ventilation, select with interchangeable elbows.

Elbow Endotracheal Tube Adapters
Ranging from 2.0/2.5/3.0/3.5/4.0mm are specially prepared for pediatric patients.

Bronchial Dilator (MDI) Port
Provides access for medication as bronchodilators into patients.

Color-Marked Tubing
Provides better recognition for insertion depth.

Model	Catheter Size	Length(mm)	Control Button	Elbow ET Tube Adapters	Type	Packaging
BT-SC1-0532NEC	05Fr/1.6mm	320	Gray	2.0mm/2.5mm	Color Open Irrigation Port	10sets/case
BT-SC1-0632NEC	06Fr/2.0mm		Light Green	2.5mm/3.0mm/3.5mm		
BT-SC1-0732NEC	07Fr/2.3mm		Purple	2.5mm/3.0mm/3.5mm		
BT-SC1-0836NEC	08Fr/2.6mm	360	Light Blue	3.0mm/3.5mm/4.0mm	Color One-way Needleless Irrigation Port	10sets/case
BT-SC2-0532NEC	05Fr/1.6mm	320	Gray	2.0mm/2.5mm		
BT-SC2-0632NEC	06Fr/2.0mm		Light Green	2.5mm/3.0mm/3.5mm		
BT-SC2-0732NEC	07Fr/2.3mm		Purple	2.5mm/3.0mm/3.5mm		
BT-SC2-0836NEC	08Fr/2.6mm	360	Light Blue	3.0mm/3.5mm/4.0mm	Color One-way Needleless Irrigation Port MDI Port	10sets/case
BT-SC5-0532NEC	05Fr/1.6mm	320	Gray	2.0mm/2.5mm		
BT-SC5-0632NEC	06Fr/2.0mm		Light Green	2.5mm/3.0mm/3.5mm		
BT-SC5-0732NEC	07Fr/2.3mm		Purple	2.5mm/3.0mm/3.5mm		
BT-SC5-0836NEC	08Fr/2.6mm	360	Light Blue	3.0mm/3.5mm/4.0mm	Color Luer Lock Irrigation Port	10sets/case
BT-SC7-0532NEC	05Fr/1.6mm	320	Gray	2.0mm/2.5mm		
BT-SC7-0632NEC	06Fr/2.0mm		Light Green	2.5mm/3.0mm/3.5mm		
BT-SC7-0732NEC	07Fr/2.3mm		Purple	2.5mm/3.0mm/3.5mm		
BT-SC7-0836NEC	08Fr/2.6mm	360	Light Blue	3.0mm/3.5mm/4.0mm	Color Luer Lock Irrigation Port MDI Port	10sets/case
BT-SC8-0532NEC	05Fr/1.6mm	320	Gray	2.0mm/2.5mm		
BT-SC8-0632NEC	06Fr/2.0mm		Light Green	2.5mm/3.0mm/3.5mm		
BT-SC8-0732NEC	07Fr/2.3mm		Purple	2.5mm/3.0mm/3.5mm		
BT-SC8-0836NEC	08Fr/2.6mm	360	Light Blue	3.0mm/3.5mm/4.0mm		

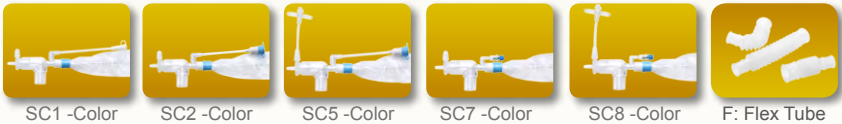

ACCESSORY CODE:(If Flex tube is desired, then include a C, F, CF in Model#.)

C: Pediatric Color catheter,
If black number coded graduations catheters are desired, then goes without "C".

BT-SCX-XXXXNECF:24 Hr Pediatric Elbow type+C(Pediatric Color catheter)+F(Flex Tube)

REMARK: Products on this page are not CE-marked.

