

# Suction Jet 72-Hour Adult T-piece

## Slide Valve Type

**BIOTEQ**® 72-Hr adult T-piece has a 90° angle for gas flow exit, it is packaged with a cap. A patient can maintain effective oxygenation delivered by the ventilator while being suctioned whether in the hospital or at home.

The catheter change frequency reduces from current every 24 hours down to 72 hours not only control the healthcare cost but also lower the cross contamination rate.

Special slide valve design seals the catheter tip into an isolated chamber, providing a thorough lavage and cleaning to prevent secretions clogging of the vacuum apparatus.

**Clear T-Piece**

Attachment to the artificial airway with 22 mm connections on both sides.

**Function switch button**

Special slide valve design seals the catheter tip into an isolated chamber, providing a thorough lavage and cleaning to prevent secretions clogging of the vacuum apparatus.

**One-way Irrigation Port**


Helps not only maintain the pressure of the ventilator circuit but also takes the injection of normal saline to clean up the catheter after each suction procedure.

**Soft Sleeve**

Highly compact and soft sleeve to effectively prevent infiltration of blood and sputum.

Model	Catheter Size	Length(mm)	Control Button	Type	Packaging
BT-SC1-1054NTD7	10 Fr/3.3mm	540	Black	T-Piece Endotracheal Open Irrigation Port	10sets/case
BT-SC1-1254NTD7	12 Fr/4.0mm		White		
BT-SC1-1454NTD7	14 Fr/4.6mm		Green		
BT-SC1-1654NTD7	16 Fr/5.3mm		Orange		
BT-SC2-1054NTD7	10 Fr/3.3mm	540	Black	T-Piece Endotracheal One-way Needleless Irrigation Port	10sets/case
BT-SC2-1254NTD7	12 Fr/4.0mm		White		
BT-SC2-1454NTD7	14 Fr/4.6mm		Green		
BT-SC2-1654NTD7	16 Fr/5.3mm		Orange		
BT-SC7-1054NTD7	10 Fr/3.3mm	540	Black	T-Piece Endotracheal Luer Lock Irrigation Port	10sets/case
BT-SC7-1254NTD7	12 Fr/4.0mm		White		
BT-SC7-1454NTD7	14 Fr/4.6mm		Green		
BT-SC7-1654NTD7	16 Fr/5.3mm		Orange		
BT-SC1-1036NTD7	10 Fr/3.3mm	360	Black	T-Piece Tracheostomy Open Irrigation Port	10sets/case
BT-SC1-1236NTD7	12 Fr/4.0mm		White		
BT-SC1-1436NTD7	14 Fr/4.6mm		Green		
BT-SC1-1636NTD7	16 Fr/5.3mm		Orange		
BT-SC2-1036NTD7	10 Fr/3.3mm	360	Black	T-Piece Tracheostomy One-way Needleless Irrigation Port	10sets/case
BT-SC2-1236NTD7	12 Fr/4.0mm		White		
BT-SC2-1436NTD7	14 Fr/4.6mm		Green		
BT-SC2-1636NTD7	16 Fr/5.3mm		Red		
BT-SC7-1036NTD7	10 Fr/3.3mm	360	Black	T-Piece Tracheostomy Luer Lock Irrigation Port	10sets/case
BT-SC7-1236NTD7	12 Fr/4.0mm		White		
BT-SC7-1436NTD7	14 Fr/4.6mm		Green		
BT-SC7-1636NTD7	16 Fr/5.3mm		Red		

ACCESSORY CODE:(If Flex tube or Disconnecting Wedge is desired, then include a W, F, WF in Model.#)  
BT-SCX-XXXXNWF72:72Hr T-Piece Slide Valve type+W(Disconnecting Wedge)+F(Flex Tube)

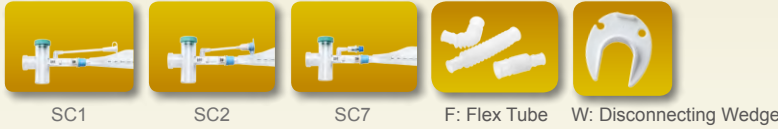

REMARK: Products on this page are not CE-marked.

