


PTCA Balloon Dilatation Catheter-Rapid Exchange

Each case is different

BIOTEQ® Rapid Exchange (Monorail) balloon meets the challenge to reach a equilibrium between various performance characteristics - pushability, trackability, torqueability, kink performance in order to provide clinical versatility in coronary artery.

When it comes to safety, everybody counts.

- Semi-compliance nylon material balloon material exhibit relatively high tensile strength, controlled compliance.
- Typical triple balloon refolding characteristics.
- The effective length of balloon catheter are available in 135cm.
- Marks on the proximal portion of the catheter shaft at 90cm and 100cm from the distal tip are for physician reference.
- Hub configuration: Enhanced strain relief design prevent kinking, transparent hub facilitates easy air bubble detection.
- Reinforced mid shaft design provides flexibility and additional support enhance overall pushability to lead physician to the target lesion.
- Guidewire compatibility: 0.014".
- Specially designed tapered tip for advanced crossability of the catheter, the tip is closely meshed with the guide wire, fish-mouth phenomenon is minimized.

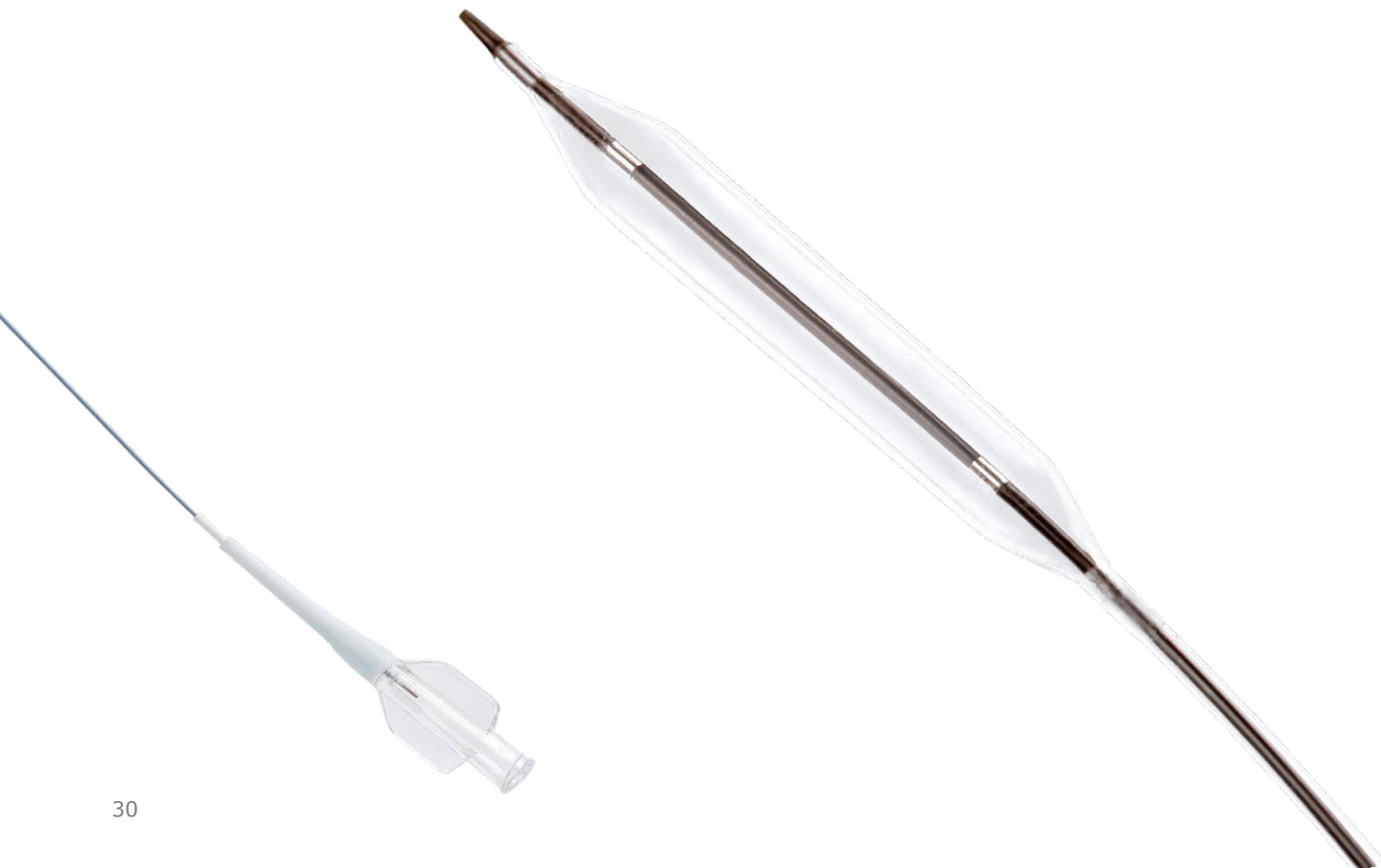

Technical Specifications

Catheter Design:	Rapid Exchange (RX)
Range of Diameter:	Ø 2.0, 2.25, 2.5, 2.75, 3.0mm
Range of Length:	60, 80mm
Guidewire Compatibility	0.014"
Usable shaft length:	135,145cm
Balloon material:	Nylon
Rated Burst Pressure(RBP):	Up to 16 ATM

Rapid Exchange Dual Lumen

Product Code	Balloon Diameter (mm)	Balloon Length (mm)	Catheter Length (cm)	Recommended Guidewire Size (inch)	Nominal Pressure (atm)	Rated Burst Pressure (atm)	Recommended Introducer Sheath Size(Fr)
BT-BC-R4-2-10-135	2.00	10	135	0.014	7	16	5
BT-BC-R4-2-15-135	2.00	15	135	0.014	7	16	5
BT-BC-R4-2-20-135	2.00	20	135	0.014	7	16	5
BT-BC-R4-225-10-135	2.25	10	135	0.014	7	16	5
BT-BC-R4-225-15-135	2.25	15	135	0.014	7	16	5
BT-BC-R4-225-20-135	2.25	20	135	0.014	7	16	5
BT-BC-R4-25-10-135	2.50	10	135	0.014	5	16	5
BT-BC-R4-25-15-135	2.50	15	135	0.014	5	16	5
BT-BC-R4-25-20-135	2.50	20	135	0.014	5	16	5
BT-BC-R4-275-10-135	2.75	10	135	0.014	5	16	5
BT-BC-R4-275-15-135	2.75	15	135	0.014	5	16	5
BT-BC-R4-275-20-135	2.75	20	135	0.014	5	16	5
BT-BC-R4-3-10-135	3.00	10	135	0.014	5	16	5
BT-BC-R4-3-15-135	3.00	15	135	0.014	5	16	5
BT-BC-R4-3-20-135	3.00	20	135	0.014	5	16	5
BT-BC-R4-2-10-145	2.00	10	145	0.014	7	16	5
BT-BC-R4-2-15-145	2.00	15	145	0.014	7	16	5
BT-BC-R4-2-20-145	2.00	20	145	0.014	7	16	5
BT-BC-R4-225-10-145	2.25	10	145	0.014	7	16	5
BT-BC-R4-225-15-145	2.25	15	145	0.014	7	16	5
BT-BC-R4-225-20-145	2.25	20	145	0.014	7	16	5
BT-BC-R4-25-10-145	2.50	10	145	0.014	5	16	5
BT-BC-R4-25-15-145	2.50	15	145	0.014	5	16	5
BT-BC-R4-25-20-145	2.50	20	145	0.014	5	16	5
BT-BC-R4-275-10-145	2.75	10	145	0.014	5	16	5
BT-BC-R4-275-15-145	2.75	15	145	0.014	5	16	5
BT-BC-R4-275-20-145	2.75	20	145	0.014	5	16	5
BT-BC-R4-3-10-145	3.00	10	145	0.014	5	16	5
BT-BC-R4-3-15-145	3.00	15	145	0.014	5	16	5
BT-BC-R4-3-20-145	3.00	20	145	0.014	5	16	5